Indications
The NICO Automated Preservation System accessories are intended for use with the Myriad for collection and preservation of resected tissue.

Device Description
The NICO Myriad System consists of a sterile, single patient use handpiece designed to be utilized only with an approved NICO console and foot pedal. Optional accessories include the Automated Preservation System, which is a collection of accessories to capture and preserve resected tissue.

The system is comprised of three parts:
 • specimen collector with filter element (collector) which collects tissue
 • specimen infusion valve (SIV) which provides the desired biological environment
 • specimen preserver (preserver) which provides the desired thermal condition

The sterile single-use Collector is attached between the Myriad handpiece and the suction canister and is used to capture tissue for evaluation. It consists of a filter element inside of an air tight housing. The filter is a removable plastic cylinder that allows the passage of fluid but captures the tissue. The collector housing is disassembled to remove the filter so that tissue may be obtained. Two styles of filters are offered, clamshell and scoop. The clamshell style opens horizontally to allow retrieval of the tissue, whereas the scoop version does not open but has a metal scoop preloaded that is removed from the filter to retrieve the tissue. Replacement or additional filter elements are offered if more than one filter is necessary during a procedure.

The sterile single-use SIV is an accessory that attaches to the Collector to allow the infusion of a desired solution directly onto the captured tissue. The SIV includes a metering line that is offered in two diameters, 0.50 mm and 0.76 mm. The metering line diameters control the rate at which fluid is delivered to the tissue.

The non-sterile single-use Preserver surrounds the Collector in an ice-chilled environment to keep the collected tissue cool. The temperature probe and color LEDs provide feedback on the condition of the ice/water bath within the Preserver.

- Refer to the relevant Myriad IFU/Manuals for additional Contraindications, Warnings, and Precautions.
- Removal of tissue using the NICO Myriad System should be performed only by persons having adequate training and familiarity with this procedure. Consult medical literature relative to techniques, complications and hazards prior to performance of any minimally invasive procedure.
- Tissue collected with the Automated Preservation System should not be re-injected into the patient.

Instructions for Use – Automated Preservation System
The components can be arranged in multiple ways depending on the desired setup. The following table, along with the figures thereafter, outline the options for setup with the Myriad handpiece:

| Myriad Handpiece | Tubing Junction on Myriad Handpiece (in sterile field) | Suction Adapter on Most Proximal End of Handpiece Tubing Prior to Waste Canister (outside sterile field) | Canister |
|------------------|--|--|----------|
| → | → SIV + Collector → | → SIV + Collector ± Preserver → | Canister |
| → | → | → SIV + Collector ± Preserver → | Canister |
| → | → SIV → | → Collector ± Preserver → | Canister |
| → | → Collector → | → Collector ± Preserver → | Canister |

- Prior to use, inspect packaging for damage or gross deformations. If there is any indication that sterility or functionality has been compromised, obtain a new device and contact NICO Customer Service about returning the original.
- Remove specimen collector and/or SIV from packaging.
- Locate desired connection point (s):
 - If in sterile field (handpiece tubing), separate blue flexible suction connector from barbed fitting and attached desired device(s); or,
 - If outside sterile field (proximal end), locate most proximal flexible suction connector and attach desired device(s).
- Connect most proximal suction connector (either from handpiece tubing or specimen collector, depending on setup) to the canister.
- If desired/applicable, prepare specimen preserver (outside sterile field) and place specimen collector within the preserver for thermal control. **Note: Ice should be added prior to assembly.** Once the lid is assembled onto the body of Preserver it is not meant to be removed. Once resection begins, pull the battery tab to activate the LEDs.
- While performing tissue resection, monitor level of tissue in specimen collector. Use additional filters/collectors as appropriate.
- Once tissue resection is complete, harvest tissue from specimen collector. Then dispose of all devices per hospital procedure.

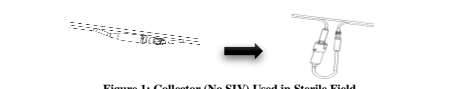
Figure 1: Collector (No SIV) Used in Sterile Field

Figure 2: Collector with SIV Used in Sterile Field (syringe not included)

Figure 3: Specimen Collector (No SIV) Used Alone Outside Sterile Field

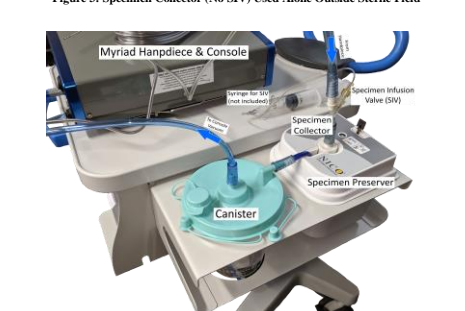
Figure 4: Specimen Collector with SIV and Preserver Used Outside Sterile Field

Operating Conditions
Temperature: +15°C to +40°C
Humidity: 30% to 75%

Storage and Handling Conditions
Temperature: -29°C to +60°C
Humidity: 30% to 80%

How Supplied
The specimen collector with filter element, replacement filter elements, and specimen infusion valves are all supplied sterile and packaged for single-use. These items are sterilized using gamma irradiation. Discard into an appropriate container after use. The specimen preserver is single-use and supplied non-sterile.

Device Related Serious Incident Reporting
Any serious incident or serious adverse event involving these products should be reported to NICO or the local Distribution Partner immediately. Serious incidents may also need to be reported to the local Regulatory Authority (e.g., Competent Authority of the Member State in EU). To report a serious incident, contact the local Representative or the Customer Service at 317-660-7118 ext.100.

Major Risks from Risk Management Activities
 • Infection
 • Transmission of Disease

Basic-UDI: 081674402MY4G2

For More Information
For more information on the NICO Myriad, contact your NICO representative or NICO Customer Service at 317-660-7118 ext.100. Additional information may be found at www.niconeuro.com, including electronic versions of the labelling: www.niconeuro.com/labelling

Patent: www.niconeuro.com/patents
© 2022 NICO Corporation. All Rights Reserved. Printed in the USA.

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© 2022 NICO Corporation. Reservados todos los derechos. Impreso en los Estados Unidos.